KÖNIG
HEALTH CARE

HC-SM10

READ THIS INSTRUCTION MANUAL BEFORE USE

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ENGLISH

Features:

Output of 2 channels -A&B can be used to treat various parts of the body simultaneously. You can choose any of the 5 treatment modes (tapping, modulation, kneading, program 1 and program 2). The wave form and frequency will fluctuate automatically in 2 computerized programs for an even more comfortable massage.

Use the HC-SM10 slimming massager for:

Stiffness of shoulders
Recovery of fatigue
Neuralgia
Paralysis of peripheral nerves
Poor blood circulation
Gradual release of muscular pain

Important preliminary information:

Warnings

Refrain from using in the following situations:

Those people who suffer from any of the following diseases or are showing any symptoms should not use the electronic pulse massager for low-frequency treatment.


1. **Heart disease**
2. **Users of pacemakers**
3. **High fever**
4. **Acute (painful) diseases**
5. **Abnormal blood pressure**
6. **Menstruation, pregnancy or after giving birth**
7. **On or around cuts and / or skin diseases**
8. **Malignant tumors**
9. **When the body is wet from bathing / showering or sweating**
10. **Tuberculosis diseases**
11. **People receiving medical treatment should consult with physician prior to using the massager.**

Caution:

1. If used by 2 persons at the same time, please ensure that each user uses only pads connected to one channel. It is advisable to always start with a low intensity setting with gradual increments to prevent any discomfort.
2. Affix the pads so that they do not overlap one another.
3. Be sure not to move the electrode pads to another part of your body or to another person's body without first turning off the power.
4. Always keep the electrode pads clean, do not touch the adhesive surface with your fingertips.
5. Should you feel any abnormality in either the main unit or your body, stop using immediately.
6. Those persons, children for example who are not capable of expressing themselves should not use the massager.
7. Do not place the electrode pads either in the mouth, on or around the genitals.
8. Do not use the unit in a room with high humidity (bathroom, etc). Or while showering/bathing.
9. Refrain from using when driving.

10. Do not use while sleeping.
11. After use, unplug the electrode cord from the main unit and coil it around the electrode holder for storage.

Body:

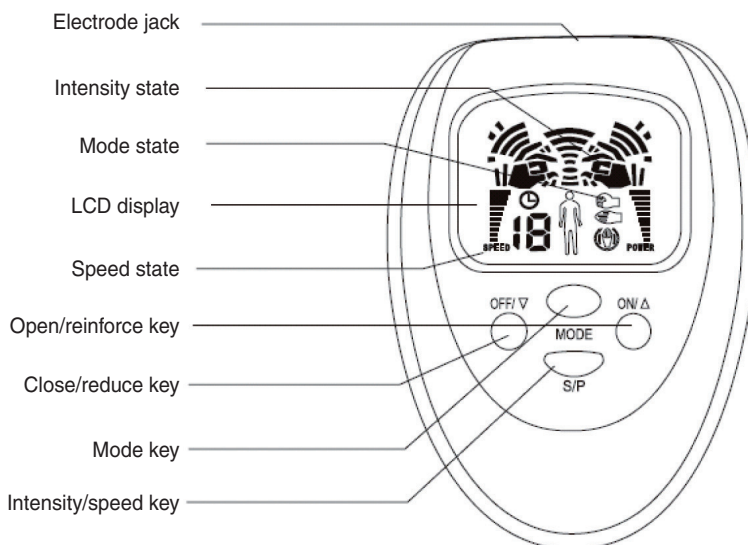Correct use:

Installing the battery

1. When replacing the batteries, be sure that the power is off.
2. Open the battery compartment by pushing it gently in the direction of the arrows.
3. Insert the battery and make sure that the (+) and (-) polarities match the diagram indicated at the back of the unit. Slide the battery cover completely shut.
4. Batt.2 x AAA (3V)

Preparing & setting the electrodes:

1. Correctly connect the electrode cord to the electrode pads.
2. Insert the plug of the electrode cord into the unit's electrode jack.
3. Peel off the film from the electrode pads.
4. Use a damp towel to wipe the area of skin on which you will be affixing the pads to remove any oil, cosmetic or dirt.
5. Should the electrode pads be soiled, both their adhesiveness, and the number of times of possible use will be reduced.
6. Apply the electrode pads to the area of the body you wish to massage. The unit cannot operate unless both pads from at least one channel are applied.

Caution:

1. Please remove batteries if the product is not in use for a long time.
2. Never apply the electrode pads to your skin with the power supply on. Doing so may result in sudden shock.
3. There are some areas of the body to which the pads cannot be easily affixed. For this case, use surgical tape from your local pharmacy to secure the pads.

Treatments:

1. Turn the power on by pressing the power switch. The TAP symbol will appear on the LCD screen.
2. Select the desired treatment mode by pressing the MODE button.
3. Raise the intensity gradually by twisting the intensity control for both channels A and B. Both channels can be set separately according to your needs.
4. When you change the treatment mode, the stimulation of the current treatment will stop. Set the intensity for the new mode by twisting the intensity control again. After treatment, turn off the power button.

Examples of treatment modes:

1. Tap - It promotes blood circulation by several kinds of stimulation such as high speed, low speed etc.
2. Modulation - It releases fatigue by several kinds of stimulation such as low speed, high speed and alternation of these speeds.
3. Kneading - It gradually releases neuralgia or muscular pain by stimulation such as strong and soft pressures.
In the above 3 modes, the stimulation in A and B channels are constant and it is possible to treat 2 part of the body at the same time.

Standard intensity, time, numbers of uses:

Intensity: The intensity of the stimulation is felt differently by each individual and also will vary with physical conditions. The strength of the batteries and the conditions of the electrodes also play a role. Therefore adjust the intensity control button to a position you feel most comfortable.

Caution:

Do not assume that by increasing the intensity of the stimulation that there will be a stronger effect. This may very well not be the case. However it is possible that excessive stimulation to the skin will result in irritation or red eruptions.

Time: A standard amount of use per area is 10-15 minutes. A single treatment should not exceed 30 minutes.

Number: A standard number of use per day is 1-2.

There may be some people who feel fatigue. In this case lower the intensity, shorten the time, and limit the number of uses to once a day.

Storage:


1. Unplug the electrode cord from the electrode jack.
2. When you have finished using the unit, set the electrode pads on to the electrode holder with the electrode cord still connected.
3. Coil the electrode cord around the groove of the electrode holder and set the plug in one of the four notches.




Maintenance and storage of the unit:

1. Clean the unit with a damp cloth or a neutral cleanser and then wipe with a dry cloth. Do not use volatile liquids such as benzene, thinner or gasoline for cleaning.
2. Do not use or store the unit where there are magnetic fields or electric wave (near TV sets or speakers).
3. Do not place the main unit in areas of high temperature, high humidity, or under direct sunlight.
4. Store the unit where there is no moisture.
5. Keep out of reach from children.
6. Remove the batteries if the unit will not be used for an extended period of time.




Maintenance and storage of electrode pads:

1. If the electrode pads become soiled the adhesive power may decrease and the skin might become irritated. If this should happen moisten the surface of the pads with water and wipe away the soiled portion, this will allow a temporary restoration of the adhesive power of the pads.
2. Do not place electrode pads in areas of high temperature, high humidity or direct sunlight. When buying replacement electrode pads, be sure to specify HC-SM10 slimming massager electrode pads for HC-SM10 slimming massager.

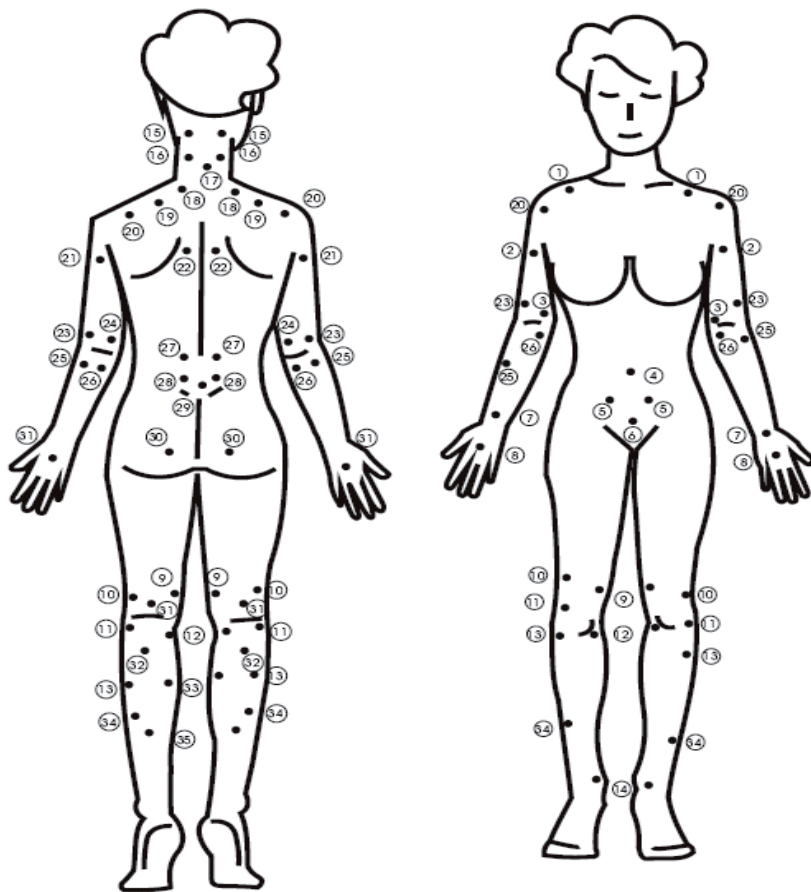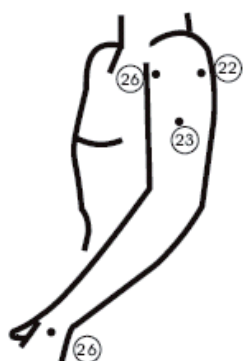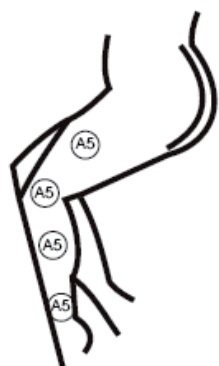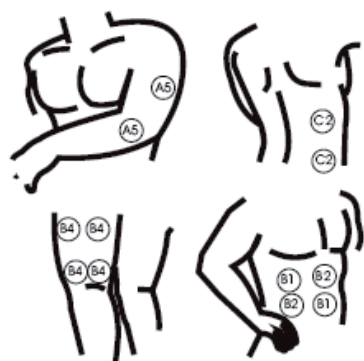
Muscle toning chart (men)

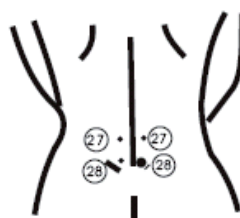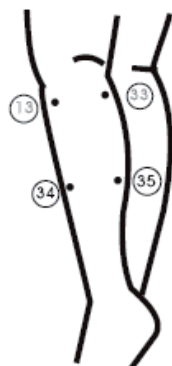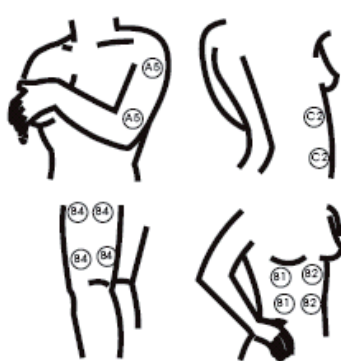
Muscle toning chart (women)


Tennis elbow

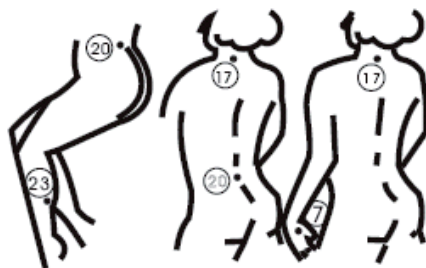
Paralysis (upper/lower body)

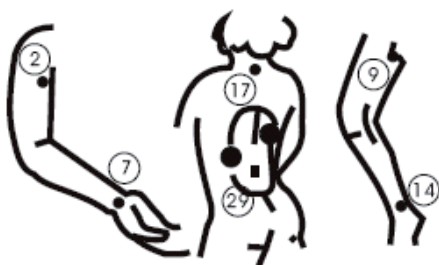
Arthritis pains on joint areas

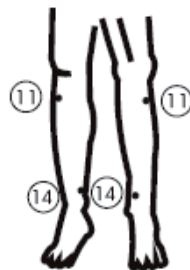
Diabetes

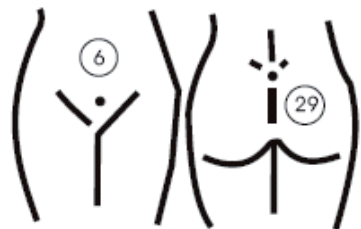
Constipation

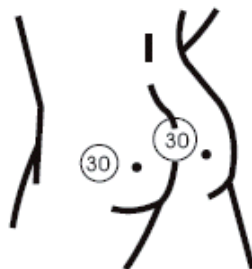vigor improvement

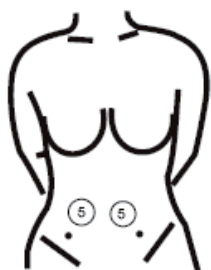Diarrhea/indigestion

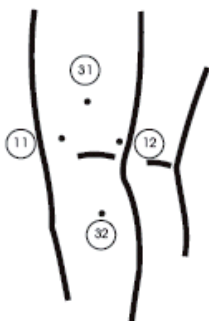
Muscles pain on limp

Knee inflammation

Rheumatoid of knee joint

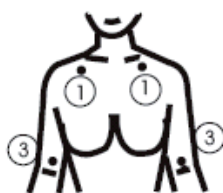Rheumatism (shoulder or elbow or any pain point)

General nervous headaches hypertension

Tinnitus

Stiff shoulder

Stiff neck

Skin turns red

1. Is the treatment time too long?
Limit your treatment time to 10-15 minutes.
2. Have the pads dried?
Wet the surface of the pads and try again. If that does not solve the problem or if the pads are quite soiled, replace them with new ones.
3. Are the pads affixed firmly?
Set them so that they fit closely to your skin.

CE DECLARATION OF CONFORMITY

We,
Nedis B.V.
De Tweeling 28
5215MC 's-Hertogenbosch
The Netherlands
Tel.: 0031 73 599 1055
Fax.: 0031 73 599 9699
Email: info@nedis.com
Internet: www.nedis.com

Declare under our responsibility that the product:

Name: König Electronic
Model: HC-SM10
Description: SLIMMING MASSAGER

Is in conformity with the following standards:

EN 55014-1:2006, EN 61000-3-2:2006
EN 61000-3-3: 1995+A1:2001+A2:2005
EN 55014-2:1997+A1:2001

Following the provisions of the EMC 2004/108/EC directive.

's-Hertogenbosch, 17-05-2010

Mrs. J. Gilad
Purchase Director

Safety precautions:

To reduce risk of electric shock, this product should ONLY be opened by an authorized technician when service is required. Disconnect the product from mains and other equipment if a problem should occur. Do not expose the product to water or moisture.

Maintenance:

Clean only with a dry cloth. Do not use cleaning solvents or abrasives.

Warranty:

No guarantee or liability can be accepted for any changes and modifications of the product or damage caused due to incorrect use of this product.

General:

Designs and specifications are subject to change without notice.

All logos brands and product names are trademarks or registered trademarks of their respective holders and are hereby recognized as such.

**Attention:**

This product is marked with this symbol. It means that used electrical and electronic products should not be mixed with general household waste. There is a separate collections system for these products.

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Copyright ©
Copyright © (Autorská práva)